User Manual
FeelNeck3
Massager with Heat

This manual should be kept on hand for easy reference.
It is your guide to a safe and efficient operation.

Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

Danger

To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning

Warning

To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Do not carry this product by its supply cord or use as a handle.
- 4) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 5) Do not drop or insert any object into any opening.
- 6) Do not use outdoors.
- 7) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 8) To disconnect, switch off all controls and remove plug from outlet.
- 9) Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons to ensure that they do not misuse the product.

10) Keep the cord away from heated surfaces. May cause serious burns. Do not use over insensitiveskinareasinthe presence of poor circulation.

11) The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

12) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

13) The appliance is only to be used with the power supply unit provided with the appliance.

14) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

15) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.

This device is for household use ONLY and must not be used for commercial or medical purposes. If you have health concerns, consult your doctor before using this massage cushion.

The massage cushion must only be used in enclosed spaces (indoor locations).



Do NOT use this product in damp rooms (e.g. when bathing or showering).

- Do NOT wash!
- Do NOT dry clean!

Name and Function of The Machine Parts

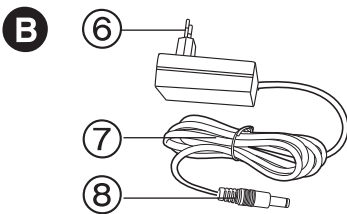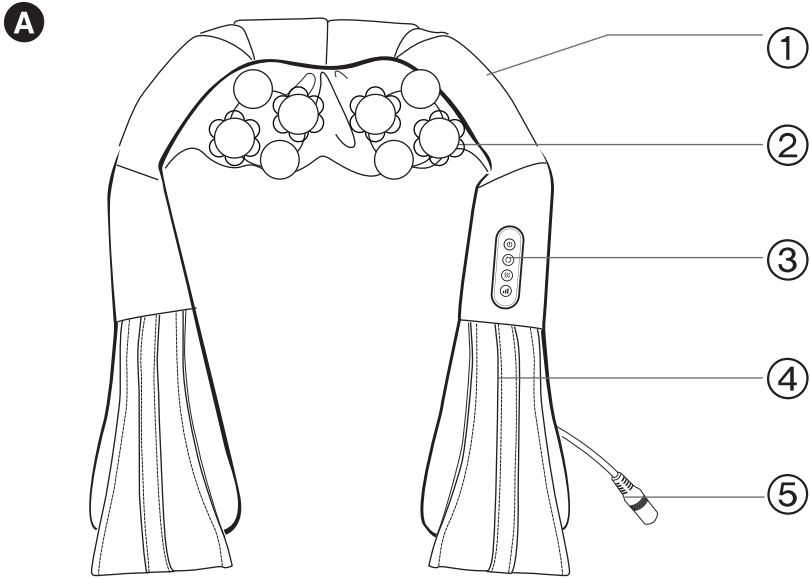A Main Unit

- ① Leather Cover
- ② Massage Head
- ③ Power Switch
- ④ Handle
- ⑤ DC Connector

B Power Adaptor


- ⑥ Plug
- ⑦ Cord
- ⑧ Appliance Plug

C Vehicle Outlet Adapter

Masssge Position

NECK MASSAGE

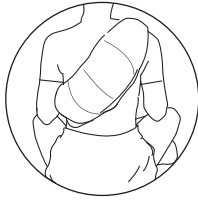BACK MASSAGE

THIGH MASSAGE


WAIST MASSAGE

ABDOMINAL MASSAGE

SOLE MASSAGE

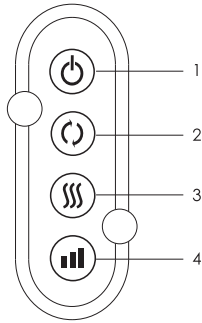How to Operate

1. Power button ON/OFF
2. Massage direction (Forward/Reverse)
3. Heating button: ON/OFF
4. Massage speed (Medium/High/Low)

How to Operate

1. Connect the power adaptor to the electrical outlet.
2. Press the power button to begin the massage. Heating and rolling function is fully activated.
3. Press the heating button to active heating function, press it again to stop.
4. Press the direction button, the massage nodes turn anticlockwise. Press it again, the direction turns clockwise.
5. Press the speed button, it will show middle / quick / slow in turn.
6. Press the power button to stop the massage at any time.

Rechargeable Description

1. Connecting the power adapter and product, insert the plug of the adapter into 100-240V ~, 50/60 Hz outlet.
 2. The massage head will appear red light flashing when it is in charging, and the light flashing will be closed when begin to working.
- ☆ During the period of using, the power adapter to charge batteries at the same time. It can automatic power off after charge completed.

Fully charge need about 5 hours

- If the battery still remains some power, the charging time will be shortened accordingly.
- After the battery fully charged, the massage head red lights flashing will be closed.

The situation of the built-in batteries has no electric power

- After massager boot, massage head can appear red light flashing, charge before use.
- Massager automatic shutdown after startup, and could not be opened, charge before use.

3. After charge completed, pull off the power adaptor from household power outlet first, then pull out the power cord plug from the power supply mouth

To Clean and Maintain the Product

1. Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket.

Please use dry cloth to clean the product.

2. Maintenance

Keep the produce clean and far from the wet environment.

Attention

Please do not use benzene, thinner and other aggressive detergents to clean the product or the product may fade.

Contact Us

If you need any assistance or support, please contact us. We will do our best to give you a viable solution.

Technology Parameter

Product: Shoulder Massager

ITEM NO. : FeelNeck3(Series JA)

Rated Voltage: DC 12V $\overline{=}$ 2.5A

Rated Power: 30W

Data of manufacture: look at massager

Made in China

	Correct Disposal of this product
	This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Авторское право © 2012г.
Ergonova GmbH